


Technical data sheet

OTICON | **Opn**

miniRITE 85
miniRITE-T 85

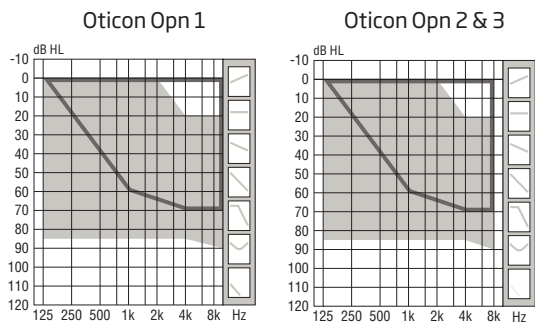

85

■ Mold, Bass & Power dome
□ Open dome



	Oticon Opn 1	Oticon Opn 2	Oticon Opn 3	
Speech Understanding	OpenSound Navigator™	Level 1	Level 2	Level 3
	- Balancing power effect	100%	50%	50%
	- Max. noise removal	9 dB	5 dB	3 dB
	Speech Guard™ LX	Level 1	Level 2	Level 3
	Spatial Sound™ LX	4 estimators	2 estimators	2 estimators
	Soft Speech Booster LX	•	•	•
Sound Quality	Speech Rescue™ LX	•	•	•
	Clear Dynamics	•	•	-
	Spatial Noise Management	•	•	-
	Fitting Bandwidth*	10 KHz	8 KHz	8 KHz
	Processing Channels	64	48	48
Listening Comfort	Bass Boost (streaming)	•	•	•
	Transient Noise Management	4 configurations	On/Off	On/Off
	Feedback shield LX	•	•	•
Personalization & Optimizing Fitting	Wind Noise Management	•	•	•
	YouMatic™ LX	3 configurations	2 configurations	1 configuration
	Fitting Bands	16	14	12
	Multiple Directionality Options	•	•	•
	Adaptation Management	•	•	•
	Oticon Firmware Updater	•	•	•
Connecting to the World	Fitting Formulas	VAC+, NAL-NL1 + 2, DSL v5.0	VAC+, NAL-NL1 + 2, DSL v5.0	VAC+, NAL-NL1 + 2, DSL v5.0
	Stereo streaming (2.4 GHz)	•	•	•
	Oticon ON App	•	•	•
	ConnectClip	•	•	•
	Remote Control 3.0	•	•	•
	TV Adapter 3.0	•	•	•
	Tinnitus SoundSupport™	•	•	•
Expected battery life, hours**	55-65	55-65	55-65	

* Bandwidth accessible for gain adjustments during fitting

** Battery size 312 - IEC PR41.

Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Oticon Opn™ miniRITE is a discreet design with a smart single push button.

Oticon Opn miniRITE-T is a discreet new style, based on the popular miniRITE, and features a telecoil and a convenient double push button for easy volume and program control.


OpenSound Navigator™ provides better speech understanding by continuously analyzing the environment, balancing all sound sources and attenuating the dominating noise.

TwinLink™ wireless technology combines binaural communication and 2.4 GHz connectivity in stereo directly to external digital devices with very low power consumption.

Oticon Opn is a Made for iPhone® hearing aid.

Oticon Opn is built on the new Velox™ platform, providing frequency resolution in as many as 64 channels (Opn 1).

Fully programmable with updatable firmware, the Velox platform is ready for the future.

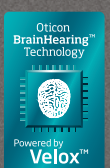

Oticon Opn is compatible with iPhone 7 Plus, iPhone 7, iPhone SE, iPhone 6s Plus, iPhone 6s, iPhone 6 Plus, iPhone 6, iPhone 5s, iPhone 5c, iPhone 5, 9.7-inch iPad Pro, 12.9-inch iPad Pro, iPad Air 2, iPad Air, iPad (4th generation), iPad mini 4, iPad mini 3, iPad mini 2, iPad mini, and iPod touch (5th and 6th generation). Devices must be running iOS 9.3 or later. Please visit www.oticon.ca/connectivity for more details on compatibility.




Technical data Measured according to		Ear Simulator IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010			ZCC Coupler ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006		
Oticon Opn miniRITE/miniRITE-T		Opn 1	Opn 2	Opn 3	Opn 1	Opn 2	Opn 3
Frequency range Hz		120-9500	120-7500	120-7500	100-8500	100-7500	100-7500
OSPL90	Peak	127 dB SPL			116 dB SPL		
	1600 Hz	120 dB SPL			111 dB SPL		
	HFA-OSPL90	121 dB SPL			112 dB SPL		
Full-on gain*	Peak	66 dB			54 dB		
	1600 Hz	52 dB			43 dB		
	HFA-FOG	55 dB			47 dB		
Reference test gain		45 dB			34 dB		
Telecoil output (1600 Hz) (miniRITE-T)	1 mA/m field	82 dB SPL			-		
	10 mA/m field	102 dB SPL			-		
	SPLITS L/R	-			94/94 dB SPL		
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 2 %			< 2 %		
	800 Hz	< 3 %			< 2 %		
	1600 Hz	< 2 %			< 2 %		
Equivalent input noise level	Omni (dB SPL)	25	26	26	20	21	21
	Dir (dB SPL)	32	33	33	29	30	30
Battery consumption**	Typical	1.6 mA			1.7 mA		
	Quiescent	1.5 mA			1.5 mA		
Battery life, artificial measurement, hours***		110			105		
IRIL (IEC 60118-13:2011) miniRITE		800/1400/2000 MHz: 31/<15/<15 dB SPL					
IRIL (IEC 60118-13:2016) miniRITE-T		700/1400/2000 MHz: 20/20/24 dB SPL					

* Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.

** Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.

*** Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994).

The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.

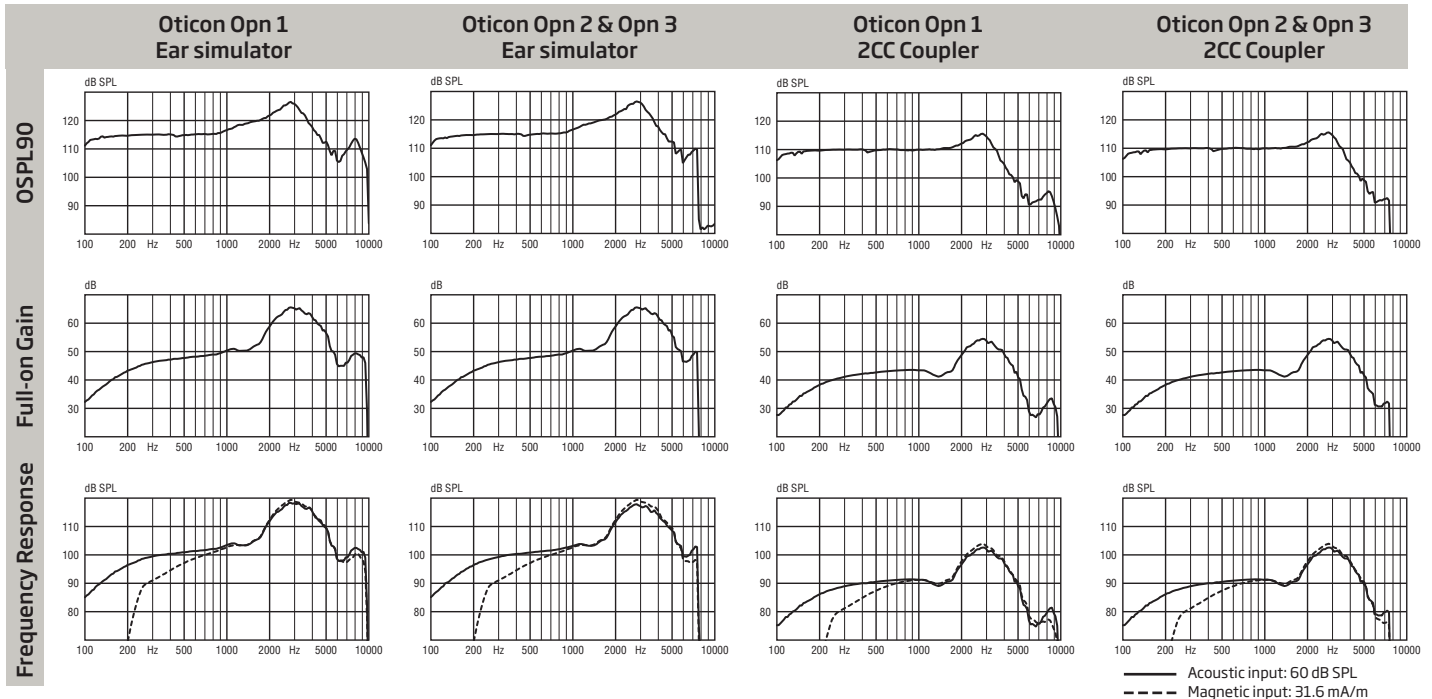

Technical information: Omnidirectional mode is used unless otherwise stated.

Operating conditions

Temperature: +1°C to +40°C

Relative humidity:

5% to 93%, non-condensing

Storage and transportation conditions

Temperature and humidity should not exceed the following limits for extended periods during transportation and storage.

Temperature: -25°C to +60°C

Relative humidity: 5% to 93%, non-condensing

Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Manufacturer:
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark
www.oticon.global

Distributor:
Oticon Canada
1-6950 Creditview Road
Mississauga, ON L5N 0A6
www.oticon.ca

oticon
PEOPLE FIRST